

GRIP-LOK[®]



Catheter and Line Securement
Safe, Simple, Secure

Table of Contents



Universal Securement..... 2

Arterial Catheter Securement 3



Foley Catheter Securement 4

Universal Securement..... 5-6



PICC Securement..... 7-8

CVC Securement..... 9-10



Specialty Securement 11-12

Hydrocolloid Securement 13



Safe... Grip-Lok is latex free and packaged sterile, allowing you to use Grip-Lok *close to the insertion site or right in the sterile field.*

Simple... Grip-Lok is applied in the same simple manner as a bandage and can be easily placed and adjusted with gloved hands.

Secure... Durably engineered, Grip-Lok holds better than tape and helps prevent dislodgement and migration, resulting in reduced risk of lost lines and patient discomfort. Using Grip-Lok for catheter stabilization conforms to the Infusion Nursing Standards of Practice¹ and can eliminate the need for sutures and unscheduled catheter restarts, thus directly reducing needlestick injuries².

Manufactured catheter stabilization devices are now listed as “preferred” in the Infusion Nursing Standards of Practice¹. Grip-Lok is a stabilization and securement device that is strong enough to lock catheters, tubes and lines securely in place, yet versatile enough for almost any securement situation.

¹Infusion Nursing Standards of Practice, February 2006, Infusion Nurses Society

²OSHA Fact Sheet: Securing Medical Catheters

Universal Securement



Universal Peripheral IV and Arterial Catheter Securement

Grip-Lok 3300MIV for winged and 3300MART for non-winged applications secure the catheter and Luer hub from the top, bottom and sides. The unique wrap-around design provides securement that maintains optimal catheter insertion angle and provides protection for the skin from hard plastic parts.

Overall length: 3.5" (8.9cm)
Patent Pending



PRODUCT #	DESCRIPTION	U/M
3300MIV	Universal Winged Peripheral IV and Arterial Catheter Securement Device, Sterile, Individually Packaged with 1 Tape Strip	100/bx
3300MART	Universal Arterial and Peripheral IV Catheter Securement Device, Sterile, Individually Packaged with 1 Tape Strip	100/bx
3300MIV-KIT-1*	IV Start Kit Sterile Packaged with 1-Grip-Lok 3300MIV, 1 Tape Strip, 1-2.37" x 2.75" Transparent Dressing, 1 Rolled Tourniquet, 1 Prevatcics 1.0mil Swab, 2 Gauze Sponges	50/bx
3300MIV-KIT-2	IV Start Kit with 1-Grip-Lok 3300MIV, 1 Tape Strip, 1-2.37" x 2.75" Transparent Dressing	100/bx


* Not available in Canada and the EU



3300MIV and 3300MART will secure ported catheters.



3300MIV and 3300MART will secure needleless connectors.

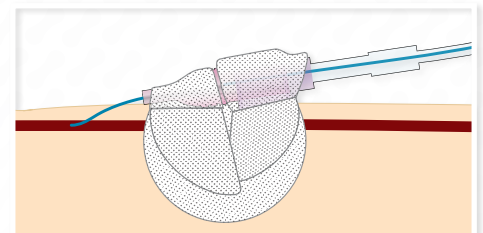

By securing both the catheter and Luer connector, the 3300MIV and 3300MART maintain insertion angle.

Arterial Catheter Securement



BD® Arterial Cannula with FloSwitch®

Grip-Lok Catheter Securement for the BD Arterial Cannula with FloSwitch has a die-cut adhesive area that matches the shape of the arterial cannula hub for proper positioning and additional horizontal securement.

Designed for BD hubs. BD and FloSwitch is a trademarks of Becton, Dickinson and Company.
Overall length: 3.5" (8.9cm)
Patent # D616,091

PRODUCT #	DESCRIPTION	U/M
3301MCS-BD*	Grip-Lok Medium for BD Arterial Cannula with FloSwitch Securement Sterile, Individually Packaged	100/bx

* Not available in the USA

Arterial Catheter


Grip-Lok for femoral and radial arterial catheters provide strong securement for introducer sheaths during cardiac catheterization procedures allowing clinicians to use both hands for more efficient manipulation of guide wires.

Overall length: 3.5" (8.9cm)
Patent Pending

PRODUCT #	DESCRIPTION	U/M
3602FRAC	Universal Femoral and Radial Catheter Securement Device, Sterile, Individually Packaged with 1 Tape Strip	100/bx



Foley Catheter Securement

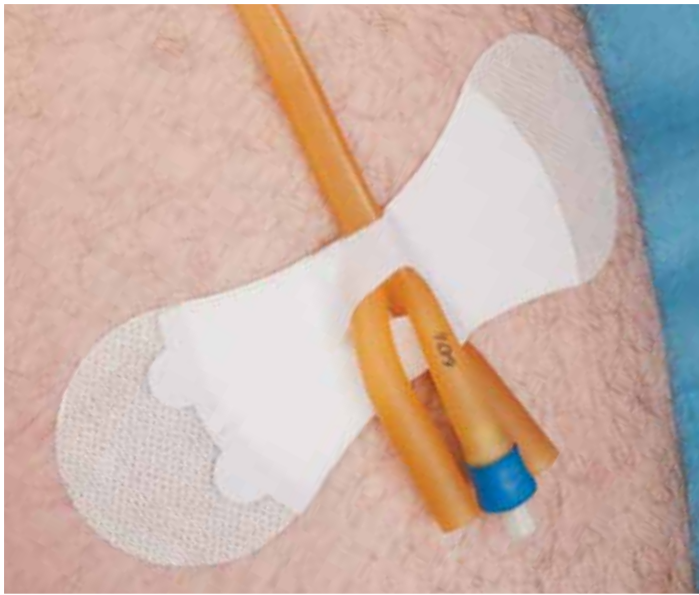

Foley Catheter

Grip-Lok Large Foley Catheter Securement is designed to hold the catheter in a channel created by a foam insert and a fabric loop over-strap. Grip-Lok for Foley Catheter does not utilize adhesive to secure the catheter in place allowing for freedom of movement and patient comfort. The degree of rotational and horizontal movement can be adjusted using the top fabric over-straps.

Designed for tubes sized 12 - 30 French in diameter.
Overall length: 6.25" (15.9cm)
Patent # D616,983

PRODUCT #	DESCRIPTION	U/M
3400LFC	Grip-Lok Large Foley Universal Securement, Sterile, Individually Packaged	100/bx


Pediatric Foley Catheter

Grip-Lok Pediatric Foley Catheter Securement works just like the larger 3400LFC, but the reduced size may be preferred for pediatric and small adult patients.

Designed for tubes sized 12 - 30 French in diameter.
Overall length: 3.5" (8.9cm)
Patent Pending

PRODUCT #	DESCRIPTION	U/M
3600PFC	Grip-Lok Pediatric Foley Catheter Securement Sterile, Individually Packaged	100/bx



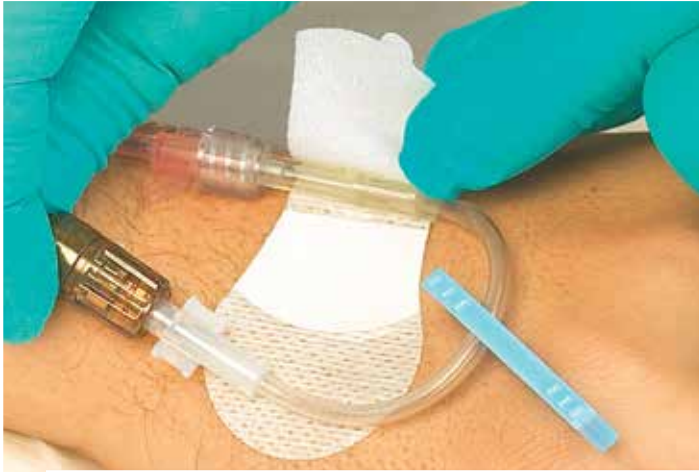

Medium

Grip-Lok Medium universal stabilization device is versatile enough for almost any securement situation and is strong enough to securely hold lines and tubes 3 to 8mm in diameter.

Designed for tubes sized 1/8" - 5/16" (3.0 - 8.0mm) in diameter.
Overall length: 3.5" (8.9cm)

Applications:

- IV Lines
- Surgical Drains
- Catheter tubes
- Bedside Monitoring Equipment
- Any tube, line or catheter sized between 3 and 8mm in diameter that needs to be secured.



PRODUCT #	DESCRIPTION	U/M
3300M	Grip-Lok Medium Universal Securement, Sterile, Individually Packaged	100/bx

Strong enough to lock catheters, tubes and lines **securely** in place, yet **versatile** enough for almost **any securement situation.**

Universal Securement

Large

Grip-Lok Large is a universal large line and catheter stabilization device that will hold tubes 6 to 13mm in diameter securely in place. Grip-Lok Large can provide firm securement for chest tubes and provides comfortable securement of wide range of latex Foley catheter sizes.

Designed for tubes sized 1/4" - 1/2" (6.0 - 13.0mm) in diameter.
Overall length: 6.25" (15.9cm)

Applications:

- Chest Tubes
- Foley Catheters
- Any tube, line or catheter sized between 6 and 13mm in diameter that needs to be secured.



PRODUCT #	DESCRIPTION	U/M
3400L	Grip-Lok Large Universal Securement Sterile, Individually Packaged	100/bx



Small

Grip-Lok Small universal stabilization device is strong enough to lock lines and tubes 1.5 to 4.5mm in diameter. Grip-Lok Small can be used in a wide variety of applications that require a smaller footprint area such as on infants or children.

Designed for tubes sized 1/16" - 3/16" (1.5 - 4.5mm) in diameter.
Overall length: 3" (7.6cm)

Applications:

- Feeding Tubes
- IV Lines
- Bedside Monitoring Equipment
- Any tube, line or catheter sized between 1.5 and 4.5mm in diameter that needs to be secured.

PRODUCT #	DESCRIPTION	U/M
3200S	Grip-Lok Small Universal Securement Sterile, Individually Packaged	100/bx



Universal PICC Securement



Medium Wide Adhesive

Grip-Lok Medium Wide Adhesive will adhere and secure catheters made from silicone, PVC and other plastic compounds. The versatility of Grip-Lok Medium Wide Adhesive will secure larger tubes in addition to a wide variety of hub sizes and shapes.

Overall length: 3.5" (8.9cm)
Patent # D608,887

PRODUCT #	DESCRIPTION	U/M
3300MWA	Grip-Lok Medium Wide Adhesive Universal Securement, Sterile, Individually Packaged	100/bx



PICC Securement

Grip-Lok CS Catheter Specific PICC Securement

“The **Grip-Lok CS series** provides a **foam die-cut insert** that **matches the shape of the catheter hub** for **proper positioning** and **additional securement stability**.”



Medcomp® Pro-PICC®



Bard™ PowerPICC™



Navilyst® Xcela® PICC



AngioDynamics® Morpheus® PICC

Designed for both single and multi-lumen catheter for the brands specified. Overall length: 3.5" (8.9cm)

PRODUCT #	DESCRIPTION	U/M
3302MCS-MC	Grip-Lok Medium for Medcomp Pro-PICC Securement Sterile, Individually Packaged	50/bx
3304MCS-BA	Grip-Lok Medium for Bard PowerPICC Securement Sterile, Individually Packaged	50/bx
3306MCS-NA	Grip-Lok Medium for Navilyst Xcela PICC Securement Sterile, Individually Packaged	50/bx
3307MCS-AD	Grip-Lok Medium for AngioDynamics Morpheus PICC Securement Sterile, Individually Packaged	50/bx

Patent # 8,500,698

Medcomp and Pro-PICC are trademarks of Medical Components Inc

Bard and PowerPICC are trademarks of CR Bard, Inc.

Navilyst, Xcela, AngioDynamics and Morpheus are trademarks of AngioDynamics, Inc.



Universal CVC Securement

Grip-Lok 3601CVC for PICC and CVC securement provides an adhesive area that will stick to a wide range of hubs and delivers superior protection from accidental line pulls.

Overall length: 3.5" (8.9cm)
Patent Pending

The 3601CVC can secure on the chest or neck and is designed to provide stable suture-free securement at the catheter attachment points.



PRODUCT #	DESCRIPTION	U/M
3601CVC	Grip-Lok Universal PICC and CVC Securement, Sterile, Individually Packaged with 1 Tape Strip	100/bx


CVC Securement



Arrow® CVC

Grip-Lok for Arrow CVC will secure subclavian access at the hub and will also fit alternative catheter attachment points.

Designed for Arrow-Howes CVC hubs. Arrow is a trademark of Teleflex Incorporated.
Overall length: 3.5" (8.9cm)
Patent # 8,500,698

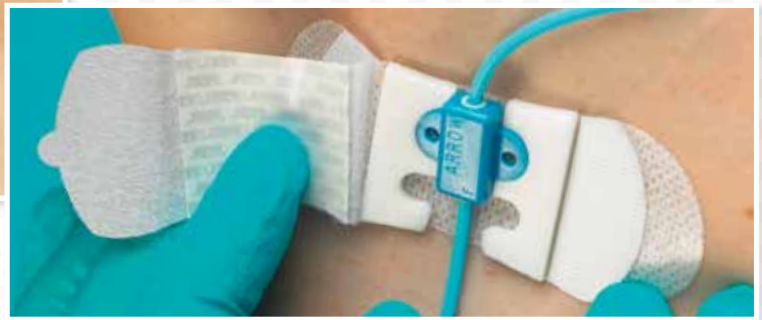

PRODUCT #	DESCRIPTION	U/M
3303MCS-TA	Grip-Lok Medium for Arrow CVC Securement, Sterile, Individually Packaged With 1 Tape Strip	50/bx

Medcomp® Dialysis CVC

Grip-Lok for Medcomp Dialysis CVC will secure the Medcomp and other popular dialysis catheters at the hub.

Designed for Medcomp Dialysis hubs. Medcomp is a trademark of Medical Components Inc.
Overall length: 3.5" (8.9cm)
Patent # 8,500,698

PRODUCT #	DESCRIPTION	U/M
3308MCS-MC	Grip-Lok Medium for Medcomp Dialysis CVC Catheter Securement, Sterile, Individually Packaged	50/bx





Nasal Gastric Securement

Grip-Lok Nasal Gastric offers firm securement of feeding tubes and will fit a wide range of tube sizes. Grip-Lok is applied in the same simple manner as a bandage and can be easily placed and adjusted with gloved hands.

Designed for tubes sized 1/8" - 13/64" (3.0 - 5.0mm) in diameter.
 Overall length: 1.75" (4.45cm)
 Patent # D663,834

PRODUCT #	DESCRIPTION	U/M
2100ANG	Grip-Lok Nasal Gastric Securement Sterile, Individually Packaged	100/bx

Large Nasal Gastric Securement

Grip-Lok Large Nasal Gastric provides a wide adhesive area to distribute the pull force across the skin surface.

Designed for tubes sized 1/16" - 1/2" (1.5 - 13.0mm) in diameter.
 Overall length: 4.25" (10.8cm)
 Patent # D720849

PRODUCT #	DESCRIPTION	U/M
3401LNG	Grip-Lok Large Nasal Gastric Securement Sterile, Individually Packaged	100/bx



Specialty Securement

Epidural Catheter Securement

Grip-Lok Epidural is a catheter stabilization device that provides a foam insert specifically designed for securing and preventing kinking of epidural catheters.

The low profile of the Medium Epidural Grip-Lok is designed for maximum patient comfort.

Designed for tubes sized 1/64" - 5/64" (0.5 - 2.0mm) in diameter.

Overall length: 3.5" (8.9cm)

Patent # 8,251,957

PRODUCT #	DESCRIPTION	U/M
3300MEPI	Grip-Lok Medium for Epidural Catheter Securement, Sterile, Individually Packaged with 1 Tape Strip	50/bx



Dialysis Fistula Needle

Grip-Lok Dialysis Fistula Needle Securement will hold the winged hemodialysis fistula needles and comes with two securements per sterile pouch.

Designed for Dialysis Fistula Needles.

Overall length: 4.5" (11.4cm)

PRODUCT #	DESCRIPTION	U/M
3500DFN	Grip-Lok Dialysis Fistula Needle Securement Sterile, Individually Packaged	100/bx*

*Two securement devices per sterile package.

Hydrocolloid Securement



Nasal Gastric Hydrocolloid

Grip-Lok Nasal Gastric Hydrocolloid provides firm securement of feeding tubes including silicone materials and will fit tube sizes ranging from 1 to 2mm in diameter. With a very small footprint, the Nasal Gastric Hydrocolloid is designed to provide securement on the smallest of patients.

Designed for tubes sized 1/25" - 5/64" (1.0 - 2.0mm) in diameter.
Overall length: 1.65" (4.2cm)
Patent # D663-834

PRODUCT #	DESCRIPTION	U/M
2100NGH	Grip-Lok Nasal Gastric Hydrocolloid Securement Sterile, Individually Packaged	50/bx




Neonatal PICC and Line Securement

Grip-Lok Neonatal PICC stabilization device will fit a wide range of hubs made from silicone and other plastic compounds. Grip-Lok is applied in the same simple manner as a bandage and can be easily placed on skin with gloved hands. 2200NUZA has zinc oxide skin contact adhesive that is gentle on sensitive skin and neonatal patients.

Designed for neonatal sized PICC lines.
Overall length: 2.25" (5.7cm)


PRODUCT #	DESCRIPTION	U/M
2200NUZA	Neonatal Universal Zinc Oxide Adhesive Neonatal PICC and Line Securement	100/bx

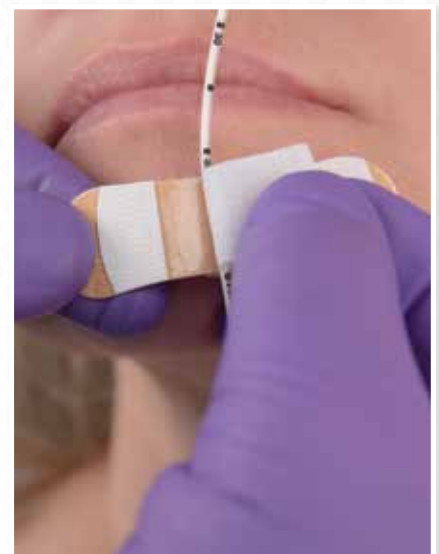


Grip-Lok products are available through your preferred distributor.


Product availability may vary by country/region. Please contact TIDI® Products for product, pricing and availability information in your area.

Visit grip-lok.com for **trial** and **inservice support** and to download **Product Brochures, Application Sheets, Instructions For Use, and Instructional Videos.**

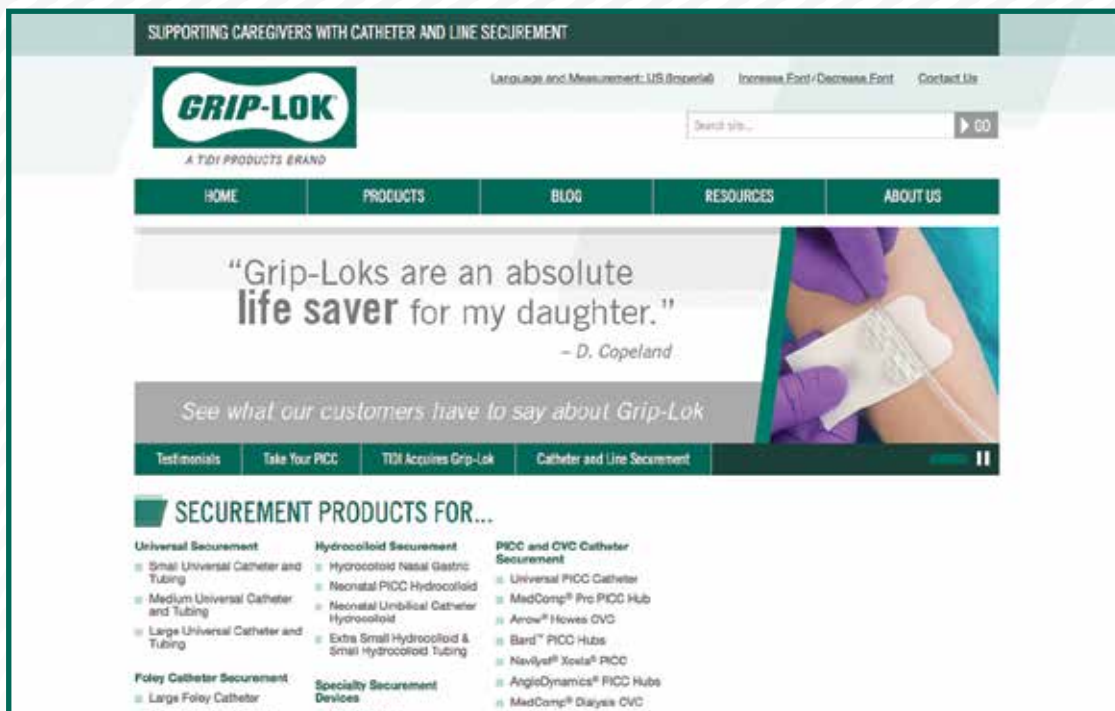

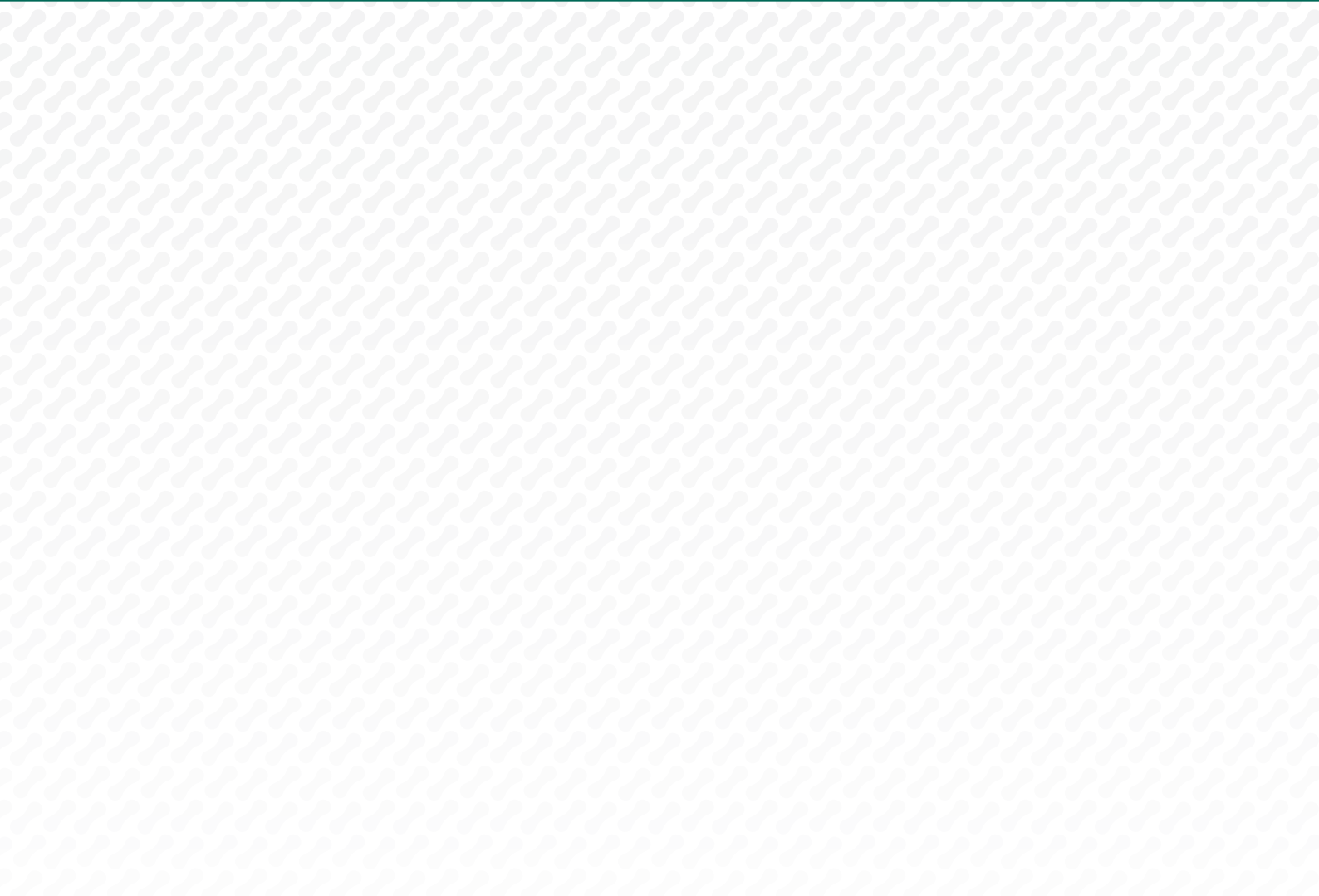

Distributed by:

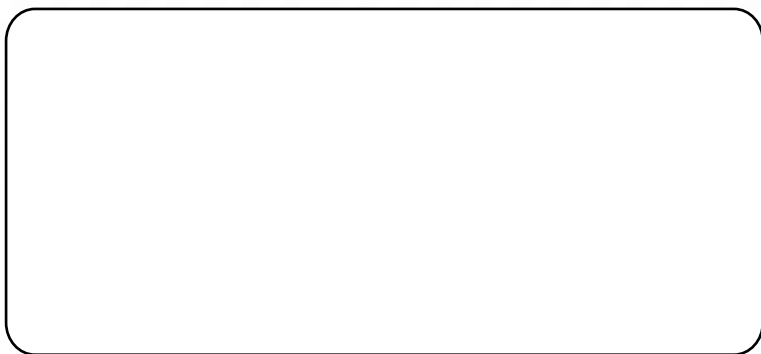

Grip-Lok is a registered trademark of TIDI Products, LLC

HYPOALLERGENIC



Not made with natural rubber latex

EC REP Authorized European Representative

Blue Box Medical Limited
Unit 29, New Forest Enterprises Centre,
Rushington Business Park, Chapel Lane
Totton, Southampton, SO40 9LA, U.K.
Tel: 023 8066 9000; Fax: 023 8066 9002



570 Enterprise Drive, Neenah, WI 54956 USA
Phone: +1 920 751 4300; Fax: +1 920 751 4370
www.grip-lok.com

MM00027 Rev.16 06/2016